


STOOL MANAGEMENT SYSTEMS

PRODUCT CATALOG



primed[®]

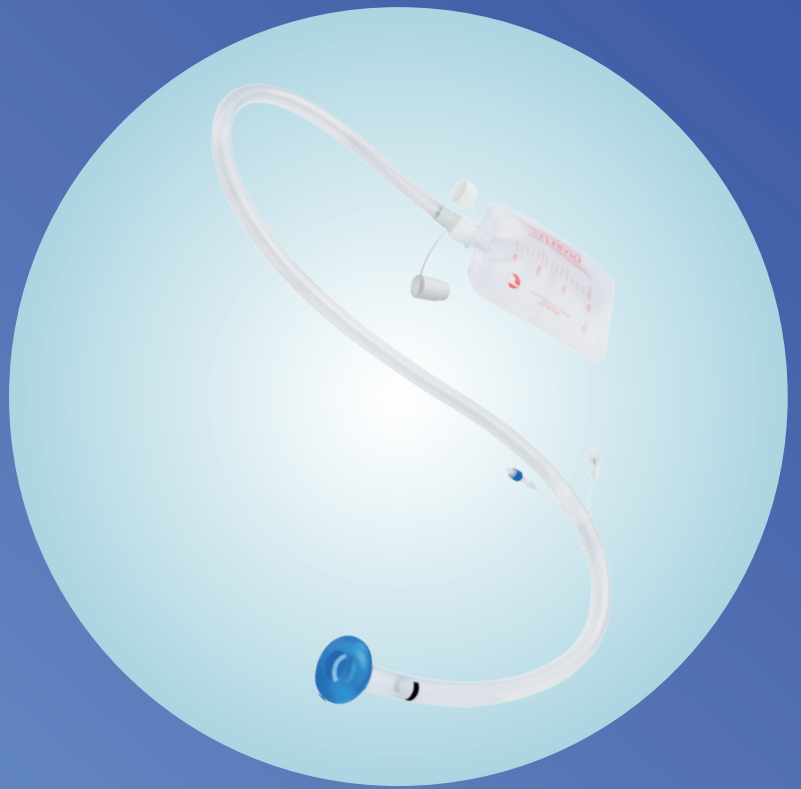

DIARFLEX®

The economic Stool Management System

Catheter hose

- Made of soft, biocompatible Silicone


Rotary closure

- For securely attaching a collection bag




Low Pressure Balloon

- Assist in holding catheter securely in place in the rectal vault
- The soft, biocompatible material prevents irritation of the sphincter muscle

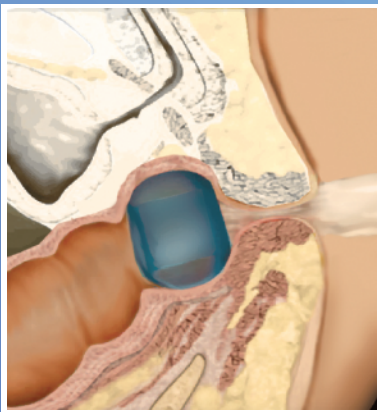

- Permits easy insertion of the catheter in the unblocked state (therefore generously apply lubricant)

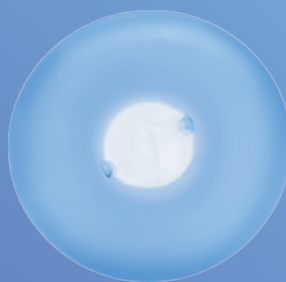

- In the blocked state, guarantees secure positioning of the catheter in the rectal opening
- Reduces the risk of leaks and liquid escaping

PRACTICAL & EFFICIENT

- Hygienic discharging and collection of liquid, semi-liquid and infectious stools
- Reduces contact with infectious microorganisms both for the nursing staff administering treatment and the patient
- Reduces dermatological complications, because infectious stool can cause severe skin damage and even lead to decubitus
- Offers the possibility to flush the device if needed
- Reduces the development of unpleasant odours
- Guarantees dignified care of the patient by eliminating the need for diapers
- Enhances the efficiency of care provision (savings on cleaning costs and working time)
- Free from LATEX
- Free from DEHP



Connector for the Low Pressure Balloon (blue, „45 ml“)

- For filling the balloon with water or saline solution
- To empty the balloon

Irrigation connector (white, „Irrig.“)

- Used to flush the device if needed

The **black marker** provides better positioning control



DIARFLEX® ADVANCED

The innovative Stool Management System

As refinement of our DIARFLEX® Stool Management System we developed DIARFLEX® ADVANCED.

Advantages:

- Offers the possibility to flush the device if needed and to take samples via the sampling port
- It is possible to rectally administer medication via the irrigation connector. For this purpose, the catheter clamp (included in delivery) enables closing the catheter system for a period of time determined by the physician
- Special odour-proof coating
- Includes hose plug for temporary seal of catheter hose during replacement of collection bag

Catheter clamp

- For closing the catheter system during rectal administration of medication
- included in the scope of delivery of DIARFLEX ADVANCED®

ADDITIONAL FEATURES



Rotary closure with hose plug

- Hose plug for temporary seal of catheter hose during replacement of collection bag



Catheter hose

- with odour-proof coating

Irrigation connector

(white, „Irrig.“)

- Used to flush the device if needed and for administration of medication



Sampling port

- For taking stool samples
- For irrigating the catheter



Low Pressure Balloon

- soft, biocompatible silicone balloon without recess between catheter hose and balloon





Collection bag for Diarflex® and Diarflex® Advanced

- To safely collect and dispose potentially infectious stools
- With superabsorbent polymer for the absorption of liquids
- Easy to connect to the hose end by means of a rotary closure
- With a firmly attached closure cap
- Capacity up to 1500 ml
- Available with an additional active carbon filter (included only in Diarflex® Advanced Stool Management Systems)

Replacement Collection Bags can be ordered separately

Additionally Included in the scope of delivery of both STOOL MANAGEMENT SYSTEMS:

- Two retaining facilities with velcro fastener for easy fixing of the system to the patient's bed
- Two replacement bags
- One syringe with Luer Lock connection for filling/ emptying the ballon and for flushing the system
- Warning labels for Fill lumen and Rinse lumen
- Instructions for use

Article number	Description	Packaging unit
20123	DIARFLEX® Stool Management System	1 set
20762	DIARFLEX® Advanced Stool Management System	1 set
20865	DIARFLEX® Replacement Bags with active carbon filter	2 pieces
20128	DIARFLEX® Replacement Bags	3 pieces

Non sterile product (usable up to 5 years)

MANUFACTURER

Primed® HALBERSTADT
 MEDIZINTECHNIK GMBH
 STRASSE DES 20. JULI 1
 D-38820 HALBERSTADT
 PHONE: 0049 (0) 3941 - 668 6
 FAX: 0049 (0) 3941 - 245 65
 WWW.PRIMED-HALBERSTADT.DE
 PRIMED@PRIMED-HALBERSTADT.DE

DISTRIBUTOR

